

# MIT Covid-19 Face Shield Folding Instructions

The Secretary of Health and Human Services (HHS) has declared a public health emergency, giving rise to [liability immunity](#) under the PREP Act for medical countermeasures against COVID-19 that satisfy certain conditions.

MIT MAKES NO REPRESENTATIONS AND EXTENDS NO WARRANTIES OF ANY KIND, EXPRESS OR IMPLIED, REGARDING THE FACE SHIELDS OR THEIR DESIGN, AND MAKES NO EXPRESS OR IMPLIED WARRANTIES OF ANY KIND, INCLUDING, BUT NOT LIMITED TO, ANY WARRANTIES OF MERCHANTABILITY, TITLE, FITNESS, ADEQUACY OR SUITABILITY FOR A PARTICULAR PURPOSE, USE OR RESULT; ANY WARRANTIES OF FREEDOM FROM INFRINGEMENT OF ANY DOMESTIC OR FOREIGN PATENTS, COPYRIGHTS, TRADE SECRETS OR OTHER PROPRIETARY RIGHTS OF ANY PARTY; AND ANY WARRANTIES THAT THE FACE SHIELDS WILL PERFORM OR COMPLY WITH FEDERAL OR STATE LAWS AND REGULATIONS REGARDING MEDICAL SUPPLIES OR MEDICAL USES.

MIT shall not be liable for any claims, demands, damages, liabilities, or injuries, including, without limitation, for property damage, bodily injury, or death, or for any indirect, special, or consequential damages (the "Claims"), arising out of or in connection with the condition, performance, ownership, use of, or inability to use, the protective face shields. To the extent the Prep Act does not apply, Recipient and/or user is exclusively responsible for the inspection, storage, transfer, disposal, and/or use of the protective face shields and any Claims resulting therefrom.


The face shield is intended to protect the eyes from splashes and is intended be worn with a mask or respirator.



1  
Peel protective film from both sides of face shield.




2  
Fold top strip over and make hard creases.

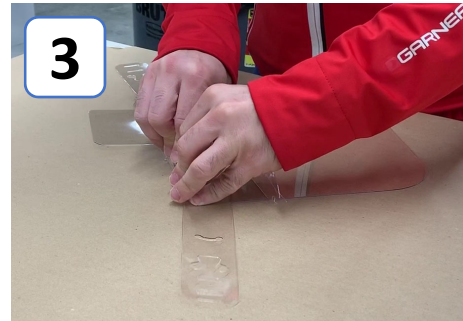

3  
Fold down visor pieces to make cover.



4  
Fold and crease side and bottom flaps

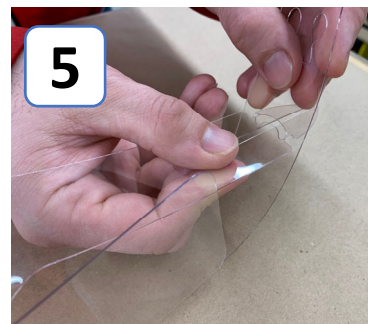

5  
Clip in tabs to the slot on side.



6  
Attach doubled hair ties to complete.



Folding video at  
<https://mit.edu/q/fs>




<https://mit.edu/q/fs>



# MIT Covid-19 Face Shield Tips & Tricks

The shield works best if the hairbands are low on the back of your head. The forehead band can go over your ears if that is comfortable for you.


You can add an extra hair tie for a looser fit or remove one for a tighter fit.

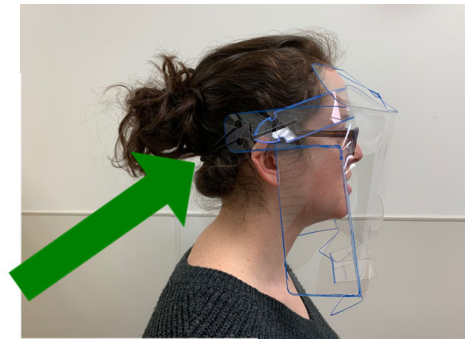

There are two resting positions of the shield--the shield should be in the down position on your face.



The up position.

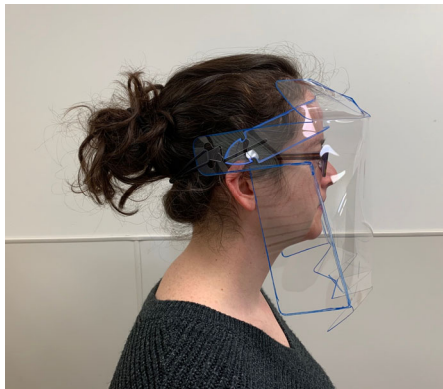

Proper down position.

For stethoscope access, feel free to open up the side edges of the shield when you pre-place your stethoscope.



The shield works best if you tightly fold the forehead piece, but loosely fold the side flaps.

